



FreeStyle

Libre 2

PLUS

Sensor Application Guide

Here, you'll discover
just how easy FreeStyle
Libre 2 Plus sensor is
to apply, activate,
and use

FreeStyle Libre 2 Plus can only be used with NFC compatible phones.

Prepare your skin

To begin, select an area on the back of your upper arm and follow these steps:



1 Wash

Use only non-moisturizing, fragrance-free soap to wash the application site and remove any oils.



2 Clean

Use an alcohol wipe to disinfect the site. If present, shave any hair on your back upper arm. Don't apply body lotions or creams until the sensor is in place.




3 Dry

Allow your skin to dry fully which maximizes stickiness. This is especially important if you plan to apply the sensor after shower or swim.

Prepare and Apply the Sensor




Peel back and remove the lid from the sensor pack.

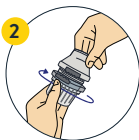

Unscrew the cap of the sensor applicator. You will hear a click.




Align the line on the sensor applicator with the line on the pack.

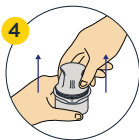

Firmly push the applicator onto the pack.



Finally, remove the applicator. The sensor is now attached and ready for application.

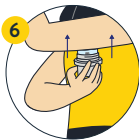

When applying, press down firmly until you hear the "click," then hold for a few seconds before pulling the applicator away



Scan to know how to use the FreeStyle Libre 2 Plus system



Important


This must be performed on a flat surface for proper sensor connection

Start the Sensor

Getting started with the LibreLink App*




Scan and download the FreeStyle LibreLink App

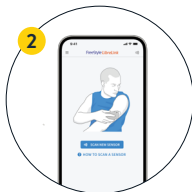

Open the FreeStyle LibreLink App* and scan the sensor. You'll receive a tone & vibration when it's activated.



Wait for 60 minutes post application.

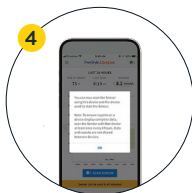

You will receive notification that the sensor is ready to go



For iOS

- **Scan readiness** : Unlock iPhone, open the App, and tap scan button in the App to put the phone in a 'scan ready' state
- **Scan technique** : Scan the sensor with the top of the iPhone
- **Report navigation** : When viewing a report, select from a drop-down menu



For Android

- **Scan readiness** : No action is required to put the phone in a 'scan ready' state when NFC is enabled
- **Scan technique** : Scan the sensor with the back of the Android phone
- **Report navigation** : When viewing a report, swipe through a 'carousel' to view other reports



Important

Make sure your Bluetooth is on, to receive data automatically every minute³

Got a question on your FreeStyle Libre 2 Plus system?
We're here to help.



1800 2100 998



adc@cs.in@abbott.com

Libre Customer Care Hours: 8:30am – 5:30pm,
Monday to Friday, excluding Saturday, Sunday and national public holidays.



Scan to visit our website
www.FreeStyleLibre.co.in

For the use of registered medical practitioners or hospitals or pharmacies
Images are for illustration purpose only. No actual patient data. Information contained herein
is for distribution outside of the USA only.

*The FreeStyle LibreLink app is only compatible with certain mobile devices and operating
systems. Please check the website for more information about device compatibility before
using the app. Use of FreeStyle LibreLink may require registration with LibreView.
Abbott Corporate Address: Abbott Healthcare Pvt. Ltd., Floor 16, Godrej BKC, Plot No. C- 68,
BKC, Near MCA Club, Bandra (E), Mumbai-400 051. Registered Address: Abbott Healthcare
Pvt. Ltd., 3, Corporate Park, Sion-Trombay Road, Chembur, Mumbai-400071.

© 2026 Abbott. All rights reserved. The sensor housing, FreeStyle, Libre, and related brand
marks are marks of Abbott. ADC-ADC-2693105 v1.0 08/26

